


Bolus

on Tandem Insulin Pumps




BOLUS

The instructions below are provided as a reference tool for caregivers who are already familiar with the use of an insulin pump and with insulin therapy in general. Not all screens are shown. For more detailed information on the operation of Tandem's insulin pumps, please visit www.tandemdiabetes.com/Support.

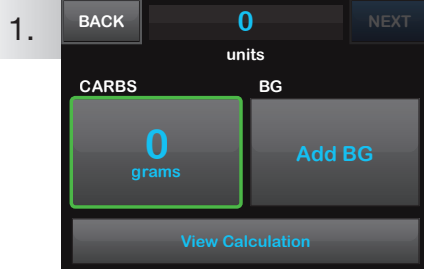


Tap **0 Grams** to enter the carbs for your bolus.

NOTE: If this button reads "units", the Carb feature is turned off in the active profile.

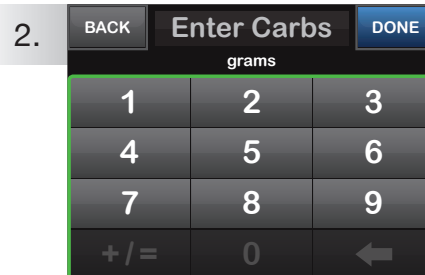


Enter desired value. Be sure "grams" is displayed above keypad for food boluses.


Tap **DONE** to continue.

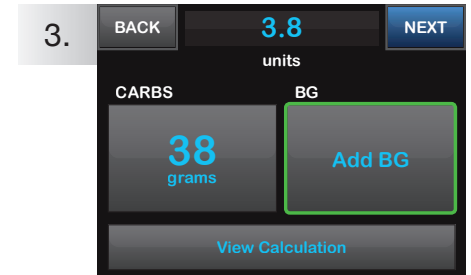


Tap **Add BG** to enter your blood glucose. Be sure "mg/dL" is displayed above keypad when entering BG values.

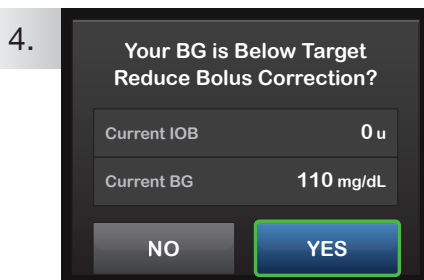


If a BG is entered that is below the target, but above 70 mg/dL, you will be offered the option to reduce the bolus amount. To accept that reduction tap **YES** otherwise tap **NO**.

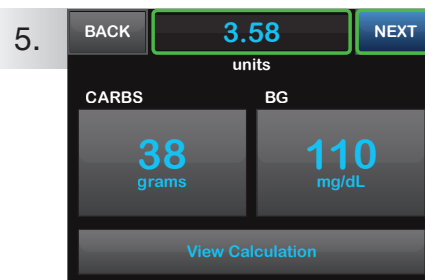


Tap **NEXT** to continue. Tap units field to manually adjust recommended dose.

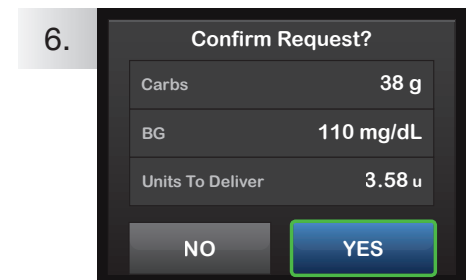


Verify the dose and tap **YES** to confirm.

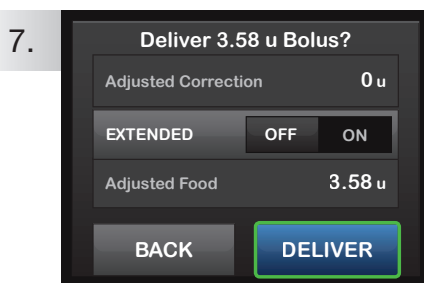

Tap **DELIVER** to deliver the food bolus immediately.

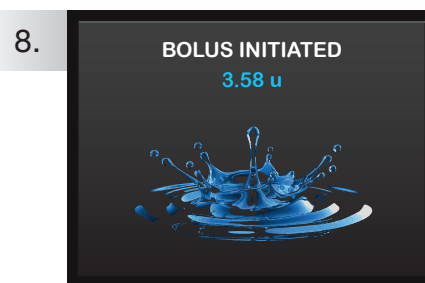

The Splash Screen confirms delivery has started.

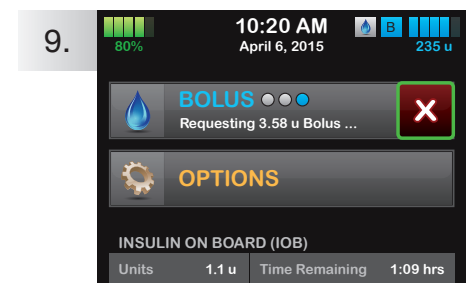

To cancel the undelivered portion of the bolus, tap the Red "X" on the BOLUS button on the Home Screen, then tap **YES** to confirm cancelled bolus.



(877) 801-6901
www.tandemdiabetes.com

A demo as simple
as our pumps!
t:simulator™ App

TRY IT TODAY!

TO DOWNLOAD APP:
Text **DEMO** to
TANDEM (826336)

