


Setting a Temp Basal Rate

on Tandem Insulin Pumps




OPTIONS

The instructions below are provided as a reference tool for caregivers who are already familiar with the use of an insulin pump and with insulin therapy in general. Not all screens are shown. For more detailed information on the operation of Tandem's insulin pumps, please visit www.tandemdiabetes.com/Support.

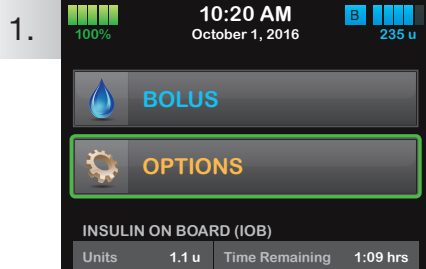


1.

Tap **OPTIONS**.

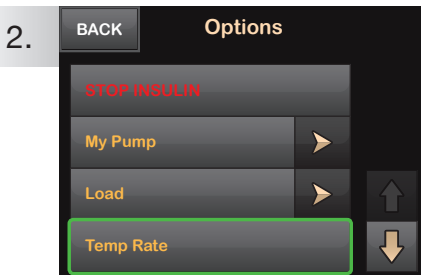


2.

Tap **Temp Rate**.

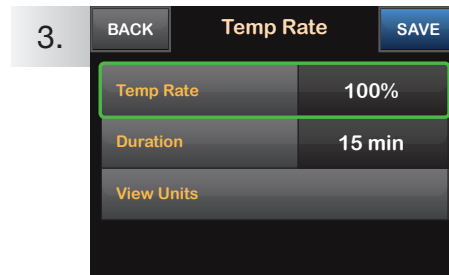

3.


Tap **Temp Rate**.

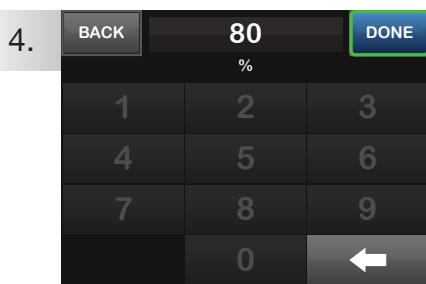

4.

Using the onscreen keypad enter desired percentage. Tap **DONE**.


NOTE: Current rate is 100%. An increase is greater than 100% and a decrease is less than 100%.

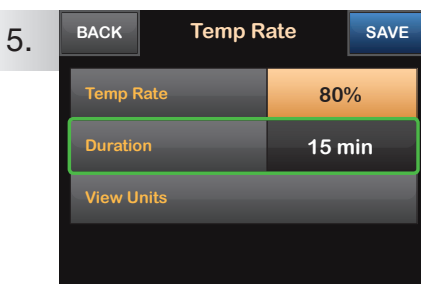

5.

Tap **Duration**. Using the onscreen keypad enter desired length of time for Temp Rate. Tap **DONE**.

NOTE: To see the actual units to be delivered tap **View Units**.

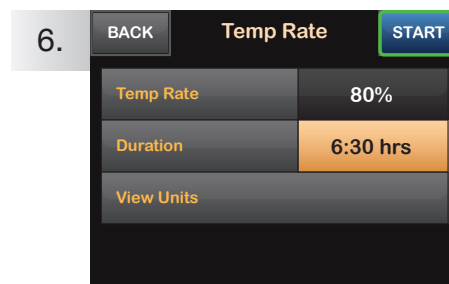

6.

Verify settings and tap **START**.



7.


The **TEMP RATE STARTED** splash screen is displayed.



8.

The Screen Lock screen will be displayed with the orange "T" icon indicating a Temp Rate is active.

NOTE: If a Temp Rate of 0% is currently active, the status of insulin delivery displays a red "T".

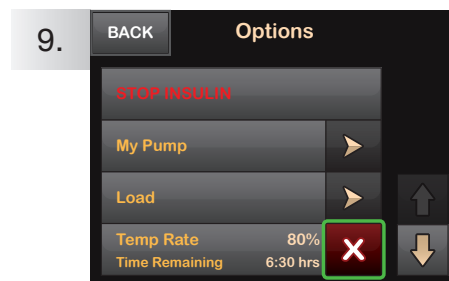

9.

NOTE: To stop Temp Rate at any time, tap **OPTIONS**, tap **X**. A confirmation screen will be displayed. Tap **STOP**.



(877) 801-6901
www.tandemdiabetes.com

A demo as simple
as our pumps!
t:simulator™ App

TRY IT TODAY!

TO DOWNLOAD APP:
Text **DEMO** to
TANDEM (826336)

